



## COMBI *touch*

### ■ stop all that stuff

There are countless instruments you currently use during an ordinary prophylaxis treatment: cures, scalers, contra angles, polishers, pastes – wouldn't it be perfect if you could have them all in one?

### ▶ start all in one

Well, now you can: the combi *touch* combines ultrasound and air polishing in one unit – it allows complete prophylaxis treatment from removal of supragingival and subgingival calculus to gentle removal of stain and biofilm and even implant cleaning.



### ► supra and subgingival


Easily switch during the treatment between supra and subgingival air polishing by simply pressing prophyl or perio button.

### ► clogging protection

Once the unit is switched on, the tubing remains clean by a consistent stream of light air.

### ► ergonomic touch panel

Due to the ergonomic touch panel, you control every function as fast and intuitive as never before – and at the same time clean and disinfect the device in literally no time.

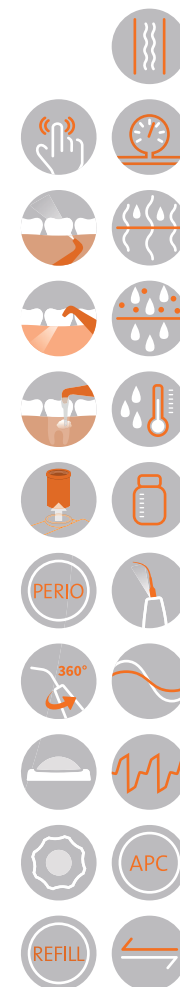

### ► led handpiece

The 360° rotatable LED cone of the ultrasound handpiece lets you set the light focus on the insert tip in every situation – regardless of the insert shape.

### ► refill function

Thanks to the exclusive „refill“ function, you can easily remove the powder chambers without switching the unit off, for fast and efficient maintenance.

### ► irrigation system

Two irrigation lines for the ultrasound system, tap water or 500 ml bottle for flexible use of water or medicinal solutions.

### ► Soft Mode

mectron's innovative SOFT MODE reduces the amplitude of the ultrasound oscillations to allow the insert to move gently and efficiently. The result: treatment efforts are met perfectly – sensitive patients are treated nearly painless and with maximum comfort.



### ► 3 different nozzles

Always best access - 90° and 120° angled spray nozzles for supra- and subgingival use in periodontal pockets up to 5 mm depth.

The Perio nozzle equipped with the disposable sterile Subgingival Perio Tip provides optimal access in pockets deeper than 5 mm.

Each nozzle is easily attached with just a click and is safely fixed to the air-polishing handpiece, thanks to a dedicated security connection system.



### ► 500 ml bottle

The illuminated bottle-system lets you use different liquids (e.g. chlorhexidine in perio or sodium hypochlorite in endo).



### ► 2 removable powder chambers

Thanks to the refill function, powder chambers are easily removable for flexible use of sodium-bicarbonate powder (prophy) and glycine powder (perio).



### ► Subgingival Perio Tip

Flexible and soft - the subgingival, disposable Perio Tip will come sterile and gets simply attached to the perio nozzle by hand.



### ► wide range of treatments

There are more than 50 ultrasound inserts available – for scaling, perio, endo and restorative treatments.



### ► Pulse Mode

The PULSE MODE transforms the ultrasound oscillation to a new profile – which results in an optimized performance for extractions and prosthetics.